

Lipsky/Schultz
Super Suction System

High Volume, Large Bore
Oropharyngeal Suction System



The LS4 is an innovative device engineered for rapid removal of upper airway obstructions and aspirated material.

Please log onto www.westmedinc.com for a complete listing of Westmed anesthesia and respiratory products.

The LS⁴™ features and benefits:

- Removes greater than 500cc of large debris within 5 seconds.
- 3/4" internal diameter for removal of large particles.
- Rigid, vent controlled, transparent wand for reliable maneuvering and accurate visualization.
- Collapse resistant corrugated tubing generates turbulence to prevent sticking.
- High Impact/Bite Resistant ABS plastic wand.
- Interchangeable Tip Guard for different suction techniques.
- Sterile Tracheal Suction Catheter Kit included with each pack.
- Quick and easy set-up and use.
- Latex Free

Westmed

The LS⁴™ Super Sucker

Lipsky/Schultz Super Suction System

High Volume, Large Bore
Oropharyngeal Suction System



Universal Tubing Adapter




Interchangeable Tip Guard and
Tracheal Suction Catheter Tip



Tip Guard



Tracheal Suction
Catheter Tip

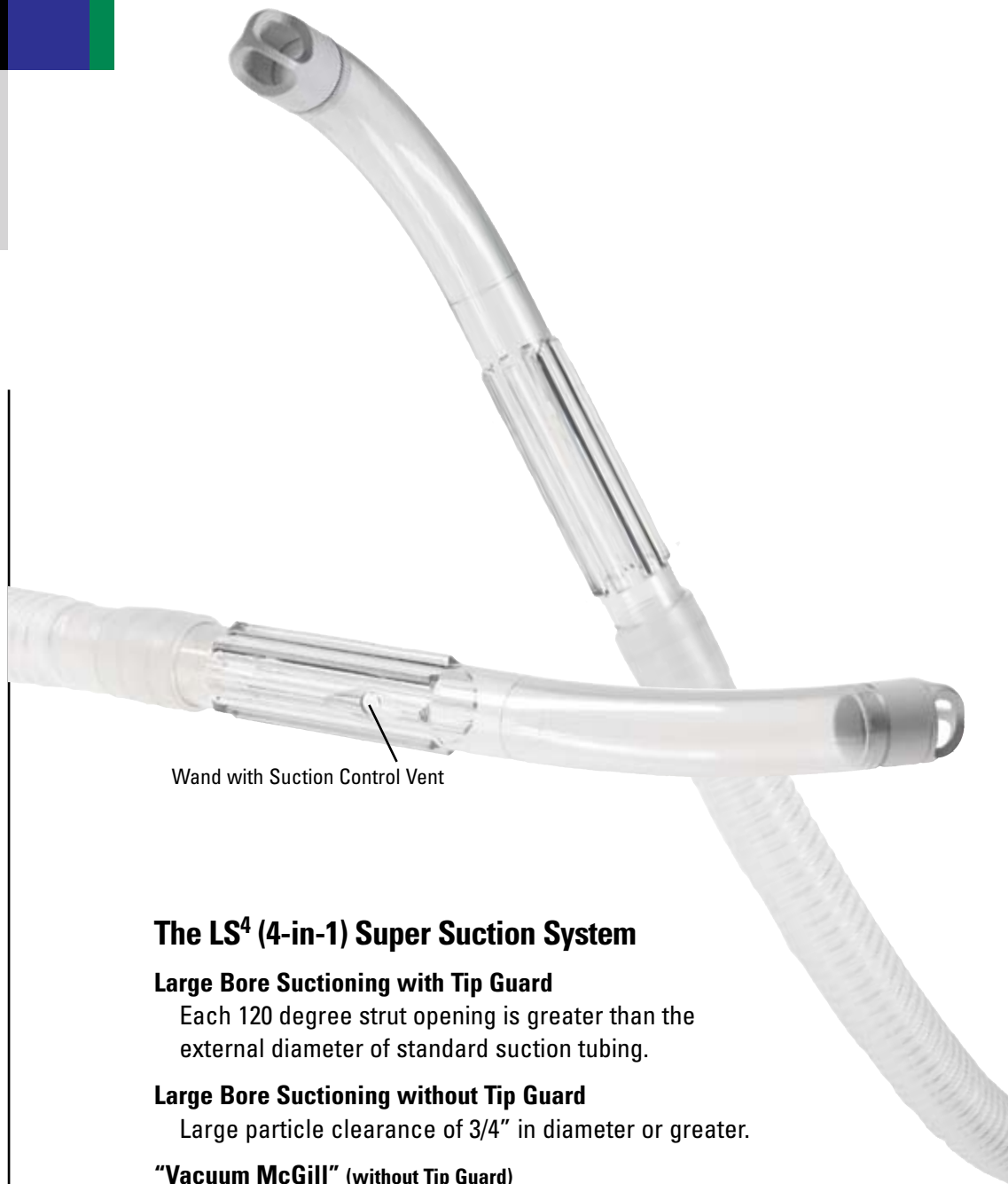

Wand with Suction Control Vent

The LS⁴ (4-in-1) Super Suction System

Large Bore Suctioning with Tip Guard

Each 120 degree strut opening is greater than the external diameter of standard suction tubing.

Large Bore Suctioning without Tip Guard

Large particle clearance of 3/4" in diameter or greater.

"Vacuum McGill" (without Tip Guard)

Can lift foreign particulate greater than 3/4" in diameter.

Flexible Sterile 16Fr Tracheal Suction Catheter Tip

For quick tracheal suctioning.

Part Number	Description	Case Quantity
Oropharyngeal Suction System		
0530	LS ⁴ , 8 Foot Tubing w/sterile tracheal suction kit	10
0531	LS ⁴ , 4 Foot Tubing w/sterile tracheal suction kit	10

Westmed